



# HaemoCer™ PLUS - Endoscopic Set

**Plant-based haemostatic powder**

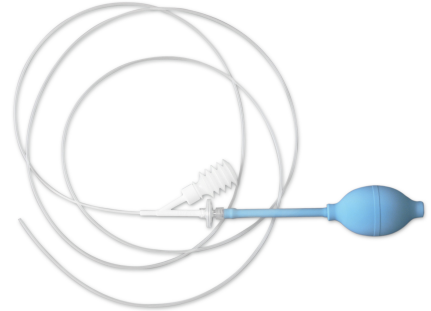

HaemoCer™ PLUS is an innovative plant-based haemostatic powder for supporting haemostasis during surgical interventions. The completely absorbable material accelerates the natural coagulation cascade due to its enormous water absorption, resulting in immediate haemostasis.

HaemoCer™ Endoscopic Applicator is used to administer HaemoCer™ PLUS in upper and lower gastrointestinal bleeding. The applicator tube is inserted through the endoscope and directed towards the bleeding. The compressed air stream generated by the pump ball transports the powder directly to the haemorrhage.

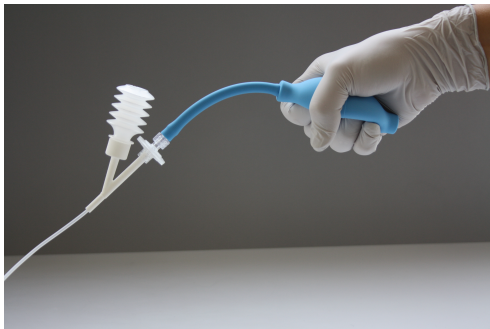

- Pump ball as a pneumatic source
- Flexible, manually operated disposable device
- **100% free** from animal and human ingredients
- **Immediately ready for use** without lengthy preparation such as thawing or mixing
- Works **without aids** such as a compressor or high-pressure gas source
- HaemoCer™ PLUS is **100% broken down** by the body through amylase **within a few days**

Article number	Description	Content per pack
HFP203	HaemoCer Plus, Haemostatic Powder, 3 g	1 unit
HFS323	HaemoCer Plus Endoscopic Set with 3 g Powder	1 unit